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
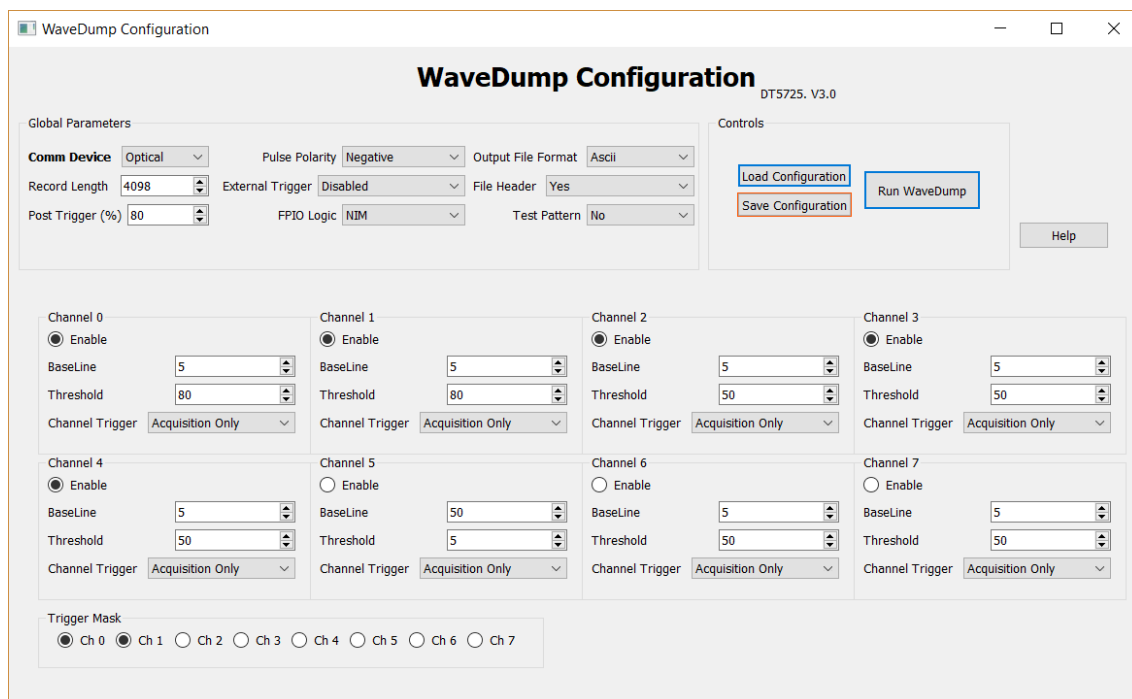
```

C:\WINDOWS\system32\cmd.exe

*****
                        Wave Dump 3.10.0 Ciemat 1.0
*****
Opening Configuration File WaveDumpConfig.txt
Connected to CAEN Digitizer Model DT5725SB
ROC FPGA Release is 04.21 - Build 3B12
AMC FPGA Release is 00.00 - Build 3A28

DAC calibration correctly loaded from board flash.

ADC Calibration check: the board is calibrated.

[s] start/stop the acquisition, [q] quit, [SPACE] help

```

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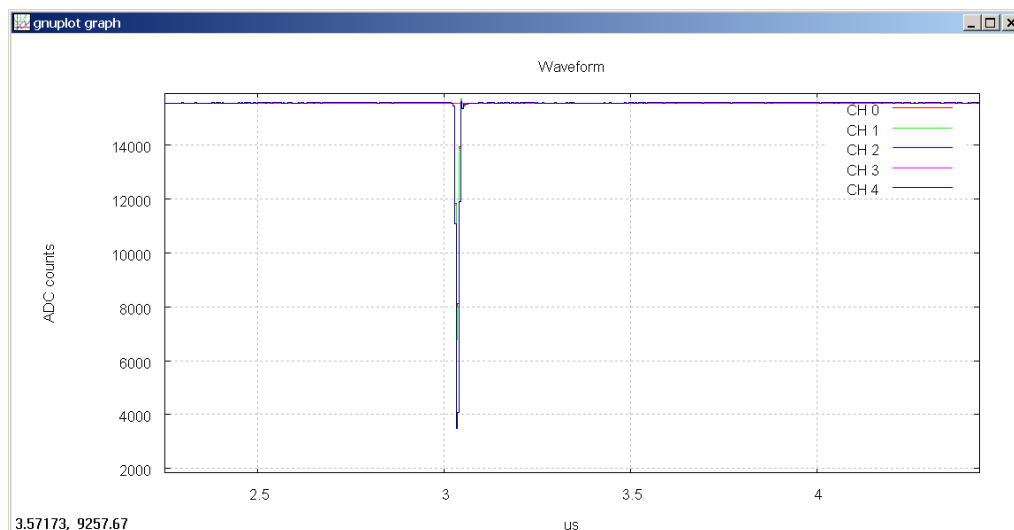Key	Function
<b>a</b>	Autoscale to be x-axis and y
<b>r</b>	Enable / Disable ruler
<b>g</b>	Enable / Disable grid
<b>y</b>	Set the scale y at full scale for the specific digitizer (scale x remains unchanged)
<b>p</b>	Return to previous zoom

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Key	Function
<b>s</b>	Start / Stop acquisition
<b>q</b>	Quit WaveDump
<b>R</b>	WaveDump restart; this command allows you to reload the configuration file and restart with a new acquisition
<b>t</b>	Single Software Trigger; this command sends a software trigger (single shot), useful especially when the card has no data (no trigger) because it forces the acquisition of an event. In analogy with the oscilloscope, this command corresponds to the “Force trigger” button.
<b>T</b>	Continuous Software Trigger; this command enables / disables the continuous generation of software trigger at a fixed rate. Inside the acquisition loop, the program sends a trigger, reads the corresponding event and executes data analysis. It corresponds to the “Auto trigger” of the oscilloscope.
<b>p</b>	Single-event plot of the waveforms of each channel enabled for plotting (which does not necessarily coincide with the n channels enabled for acquisition). The plot is single-shot.
<b>P</b>	Continuous Plot; this command enables / disables the continuous plot of the enabled channels. When enabled, the input signal trace of channel 0 is plotted by default. Check the “c” command to enable the other channels.
<b>w</b>	Save Single Event to Disk; this command causes the writing of a file for each enabled channel, named wave_n.txt, where n is the channel number. <b>Warning:</b> the files are overwritten each time. Each file can be ASCII or binary and may contain a header or not, depending on the settings assigned in the configuration file.
<b>W</b>	Continuous Event Saving; this command enables / disables the continuous events saving to file. As for the “w” command, it creates one file per channel; then it writes the events consecutively. <b>Attention:</b> in this mode, the file size can grow up very quickly.
<b>0 .. 7</b>	This command insert / remove channel n (n = 0 .. 7) from the plot (if such a channel is enabled). In case of x740, n refers to the channel within the currently active group, for example, if it's on group 2 (channels 16 to 23), press 3 to insert / remove the channel 19. In case of V1725 and V1730, n refers to channel n or n+8, depending on the visualization mode enabled by the command ‘g’ (see below).
<b>g</b>	In case of 740 and 742 series, this command switches to the next group of 8 channels (only for). <b>Note:</b> the active group refers only to the plot (in fact, the plotter can handle only 8-traces), while the acquisition is always enabled on all groups. In case of V1725 and V1730, this command switches the plot visualization from 0-7 channels to 8-15 channels and vice versa (see also the “0..7” description).
<b>f</b>	This command toggles between waveform plot and FFT plot.
<b>h</b>	This command toggles between waveform plot and histogram samples amplitude plot (not supported by 742 series).
<b>m</b>	This command displays the temperature values (in °C) of the ADC channels. It is supported only by 725, 730 and 751 digitizer families.
<b>c</b>	This command performs the channel calibration required by 725, 730, 751, and 761 digitizer families (meaningless in case of x725S and x730S digitizers).
<b>D</b>	This command performs the manual DAC calibration on the input channels (not supported by 742 series). Calibration is performed on all channels at once; the calibration coefficients are then saved into the on-board FLASH, and automatically reloaded by the software at next start-ups. Calibration is applied only if the <code>BASELINE_LEVEL</code> is used in the configuration file (see Chap. 4).
<b>[Space]</b>	This command displays the online help

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